



Excellence is Our Specialty

**MANUFACTURER & EXPORTER  
MEDICAL CONSUMABLES**



**PRODUCT CATALOGUE**

Matu International Pvt. Ltd.



# About us

**Matu International Pvt. Ltd.** is a leading manufacturer and exporter of medical, surgical and diagnostic products. Our aim, is to provide high quality health care products at affordable price to our valuable customers. We supply our product to leading and distinguished Govt. institutions, Private Corporate Hospitals, Nursing Homes, Clinics and Retail counters in India and to several other countries as well. Seamless services to our customers since a long time has helped us earn this goodwill. We thrive in complexities and work towards fulfilling the need of the hour and our customers. One of our company's main objectives is to discover and develop new and innovative products to enhance the overall health of the citizens, leaving a positive impact.

# Contents

Product Description

Page No.



## INFUSION THERAPY™

Blood Administration Set

Disposable Syringes

Infusion Set

IV Cannula

Measured Volume Set

Scalp Vein Set

Three Way Stop Cock

Three Way Stop Cock With  
Extension Tube

IV Flow Regulator

Pressure Monitoring Line



## ANAESTHESIA & RESPIRATORY CARE

Nebulizer Mask

Oxygen Mask

Twine Bore Nasal Oxygen Set

Nasal Oxygen

Suction Catheter

Spinal Needle

Endotracheal Tube

Infant Mucus Extractor

Guedel Airway

PVC Manual Resuscitator (Ambu Bag)



## UROLOGY

Folleys Ballon Catheter (Two Way)

Product Description

Page No.

3 Way Folley Catheter  
Urine Collection Bags with  
Measured Volume Set



## **Surgery and Wound Dranaige**

Pop Bandage B.P.  
Elastic Adhesive Bandage B.P.  
Paraffine Gauze Dressing  
Surgical Paper Tapes  
Sutures  
Combine Dressing  
Surgical blades  
Cotton Crepe Roll  
Yankaure Suction Set  
Closed Wound Suction Set  
Underwater Seul Drainage  
System



## **Gastroenterology**

Ryle Tube  
Feeding Tube



## **Gloves**

Latex Examination Gloves Non-Sterile  
Nitrile Examination Gloves Powder Free Non-Sterile  
Latex Examination Gloves Sterile



## **Others**

Non-Woven Surgeons Gown  
Non-Woven Doctors/Nurse  
Caps  
Cannula Fixator  
Umbilical Cord Clamp



# INFUSION THERAPY

# Blood Administration Set

Ref. no- BA0136



- Clear, Transparent & Flexible Single/Double drip Chamber
- Super Smooth kink resistant tube.
- Very efficient roller controller for accurate flow rate.
- Self-sealing Plain Latex / Bulb Latex for easy flushing.
- Needle size: 18
- EO Sterile.
- Individual Packed.

**AVAILABLE IN:** HM-Poly Pouch/Ribbon Poly Pouch/Paper Pouch/Blister Pack.

## TECHNICAL SPECIFICATION

Type	Capacity	Packing Inner/Outer
Vented	20 Drop/ml	25/500
Non-Vented	20 Drop/ml	25/500

# Disposable Syringes



- With Disposable Hypodermic Needle.
- Made from Medical grade Transparent Polypropylene.
- Sufficiently transparent barrel allows easy measurement of the volume contained in the disposable syringe friction free.
- Friction free movement of plunger ensures proper vacuum and no dead space.
- Consist of: Barrel, Plunger, and Standard Piston.
- In superior quality by high accurate moulds.
- EO Sterile

**AVAILABLE IN:** Ribbon & Blister Pack.

## TECHNICAL SPECIFICATION

Size	Ref. no.	Size	Ref. no.
2 ml	DSY01432	10 ml	DSY014310
3 ml	DSY01433	20 ml	DSY014320
5 ml	DSY01435	50 ml	DSY014350

# Infusion Set



- Clear, Transparent & Flexible Moulded Chamber (Vented or Non-Vented).
- Super smooth kink resistant tube.
- Very efficient roller controller for accurate flow rate.
- Self-sealing Plain Latex / Bulb Latex for easy flushing.
- Needle size G: 21/22.
- Available with Y-Connector & Luer- Lock.
- EO Sterile.
- Individual Packed.

**AVAILABLE IN:** HM-Poly Pouch/Ribbon Poly Pouch/Paper Pouch/Blister Pack.

## TECHNICAL SPECIFICATION

Type	Ref. no.
Economy	IV0126T1
Vented	IV0135A1
Non-Vented	IV0135A
Y-Connector	IV0135A
Luer Lock	IV0135A3



# IV Cannula



- Kink resistant FEP/PTFE Teflon with RO Cannula has extremely smooth inner and outer surface to ensure minimum trauma and smooth flow.
- Molded cannula tip ensures smooth & less traumatic cannulation.
- Injection port with non-return valve is provided for additional medication.
- Transparent flash back chamber allows immediate detection of blood.
- Color coded by size for easy identification.
- EO Sterile.
- Individual Blister packed in Medical Grade Paper & Printed Box Packing.

**AVAILABLE IN:** 14G-26G sizes

## TECHNICAL SPECIFICATION

Size	Flow Rate(water) ml/min	Box Qty.(Inner/Outer)	Ref. No.
14g	270	100/1000	IVC013814
16g	200	100/1000	IVC013816
18g	90	100/1000	IVC013818
20g	60	100/1000	IVC013820
22g	36	100/1000	IVC013822
24g	23	100/1000	IVC013824
26g	15	100/1000	IVC013826

# Measured Volume Set



- Flexible, Transparent & soft cylindrical measured volume chamber of 110ml/150 ml.
- Efficient roller controller for better flow control.
- Special Beveled strong spike for smooth penetration.
- Specially designed floating auto shut off valve acts as floating indicator and automatically shut off the drain path when the chamber gets empty to prevent air trapping in the fluid line.
- Available with or without 'Y' connection.
- Needle size G: 23
- EO Sterile.
- Soft and Kink Resistant tube. Bulb Latex or Tube latex for easy flushing.
- Separate plug for extra medication and continuous change over.
- Individual Packed.

**AVAILABLE IN:** HM-Poly Pouch/Ribbon Poly Pouch.

## TECHNICAL SPECIFICATION

Size	Capacity	Packing Iner/ Outer	Ref. No.
110ml	Dispo M.V.- 110ml	25/200	MV0137110
150ml	Dispo M.V.- 150ml	25/200	MV0137115

# Scalp Vein Set



- Designed to provide rapid venous access, for patient comfort during infusion.
- Siliconized and thin walled high quality needle for a traumatic insertion.
- Butterfly wings facilitate for better and easy fixation.
- Super smooth kink resistant tube.
- Color coded butterfly wings by size for easy identification.
- EO Sterile.
- Individual Blister or Poly Packed.

**NEEDLE LENGTH IN ALL SIZES:** 19mm

## TECHNICAL SPECIFICATION

Size	Needle Diameter	Packing Inner/Outer	Ref. No.
18g	1.20mm	100/1000	SV014118
19g	1.10mm	100/1000	SV014119
20g	0.90mm	100/1000	SV0141120
21g	0.80mm	100/1000	SV0141121
22g	0.70mm	100/1000	SV0141122
23g	0.60mm	100/1000	SV0141123
24g	0.50mm	100/1000	SV0141124

# Three Way Stop Cock



Ref. no- 3WS0118

- Three way stopcock is suitable for use during pressure infusion and invasive blood pressure monitoring.
- Clear body for easy visual of fluid flow.
- Arrow mark on tap for indicating flow direction.
- Prevent vortex and reduces thrombus formation.
- Excellent pressure resistance upto 300kpa ( 45psi).

**AVAILABLE :** Red and blue pegs for arterial and venous line identification.

# Three Way Stop Cock with Extension Tube



- Male kink resistant PVC Tube.
- Male luer lock connector at one end and female luer lock connector at other end.
- With or Without “Y” Connection.

**AVAILABLE:** with or without clamp.

**Ext.-int. Diameter :** 3.8/2.8 mm

## TECHNICAL SPECIFICATION

Size	Length(cm)	Pcaking	Ref. no.
10cm	10	50/500	3WS013910
50cm	50	50/500	3WS013950
100cm	100	30/300	3WS0139100
150cm	150	30/300	3WS0139150
200cm	200	30/300	3WS0139200

# IV Flow Regulator



- The flow regulator is used to regulate the flow of I.V. fluid from the infusion set to the intravenous cannula and having a smooth kink resistance tube, ensures a constant flow rate.
- Precise control of flow rate of IV fluids with a range of 5 to 250 ml/hr.
- With "Y" injection port (latex) for intermittent medication.
- Soft and kink resistant PVC tubing.
- Standard tube length: without "Y" site: 20cm & with "Y" site: 40cm
- Tube dia. 10 3.0 mm, OØ 4.1 mm

## TECHNICAL SPECIFICATION

Size	Ref. no.
20cm	IVFR014220
30cm	IVFR014230
40cm	IVFR014240

# Pressure Monitoring Line



- It has small bore for high pressure (Upto 140psi).
- Made of P.V.C.
- Male luer at one end female luer at other end.
- Packed in peel open pouch or soft blister pack.
- Sterile/Non-Pyrogenic/Single use.

**TUBE DIAMETER :** Internal 1.0mm & external 2.0mm

**CAPACITY :** 1.0ml per 100cm

## TECHNICAL SPECIFICATION

Size	Ref. no.
10Cm	PM014010
50Cm	PM014050
100Cm	PM0140100
150Cm	PM0140150
200Cm	PM0140200



# **ANAESTHESIA & RESPIRATORY CARE**



# Nebulizer Mask



- Designed to be used with any standard Nebulizer kit.
- Made from a soft clear plastic and are anatomically formed for comfort.
- Latex-Free.
- It has an Elastic Strap for added comfort.
- Each Mask is individually sealed and wrapped for safety.
- The Nebulizer Kit just plugs into the standard fitting located at the bottom of the Mask.
- Masks are designed to deliver more medication directly to the patient due to the direct flow design.
- Masks are contoured to fit the face of adult and are very comfortable to wear.
- Tube Length: 200 cm.
- EO Sterile

**AVAILABLE IN:** HM Poly Pack/Ribbon Poly Pack.

## TECHNICAL SPECIFICATION

Size	Tube Length	Packing Inner/Outer	Ref. no.
Adult	200cm	120	NM0110A
Pediatric	200cm	120	NM0110P
Neonate	200cm	120	NM0110N

# Oxygen Mask



- Made of clear, non-toxic PVC.
- Moulded face mask has adjustable elastic strap and integrated nose clip for proper positioning of mask on the mouth and nasal area.
- Specially designed for convenient oxygen therapy.
- Provided with long star lumen tube to ensure continuous flow of Oxygen.
- The tube has funneled shaped connector at the proximal end to facilitate simple & safe connection with the oxygen Tube.
- EO Sterile.

**AVAILABLE IN:** HM Poly Pack/Ribbon Poly Pack.

## TECHNICAL SPECIFICATION

Size	Tube Length	Packing Inner/Outer	Ref. no.
Adult	200cm	120	OM0109A
Child	200cm	120	OM0109C
Neonate	200cm	120	OM0109N

# Nasal Oxygen Cannula



- A nasal cannula is a medical device to provide supplemental oxygen therapy to people who have lower oxygen levels.
- There are two types of nasal cannulas: low flow and high flow.
- The device has two prongs and sits below the nose.
- The two prongs deliver oxygen directly into your nostrils.
- Made from soft medical grade PVC.
- Designed for easy administration of oxygen and comfort of patient.
- Star lumen main tube to avoid accidental nasal tips for maximum patient comfort.
- Twin prong nasal tips to ensure equal volume of oxygen to both air passages.
- Non-Toxic, Pyrogen free.

## TECHNICAL SPECIFICATION

Size	Ref. no.
Adult	NO0111A
Child	NO0111C
Neonate	NO0111N

# Suction Catheter



- Suction Catheter is a long flexible tube which is used in intensive case remove secretion from the mouth and airway of critically patients. This Catheter to a collection container and device that generate suction
- Manufactured from Non-Toxic, DEHP Free P.V.C. tbe.
- Universal colour coded connector for size identification and secure connection with suction connector.
- Frost finished or clear tube along with atraumatic eyes for smooth intubation.
- Graduated marked for correct position of catheter.
- Straight peel open pack.

**LENGTH :** 50cm

## TECHNICAL SPECIFICATION

Size	Ref. no.	Size	Ref. no.
6FG	SC01156	14FG	SC011514
8FG	SC01158	16FG	SC011516
10FG	SC011510	18FG	SC011518
12FG	SC011512	20FG	SC011520

# Spinal Needle



- Suitable for Spinal anaesthesia, high flow rate enables faster cerebro spinal fluid flash back.
- Quincke point for smooth penetration.
- Color coded hub for easy identification.
- Sterilization by Ethylene Oxide.

## TECHNICAL SPECIFICATION

Size	Ref. no.	Size	Ref. no.
17	SN012417	22	SN012422
18	SN012418	23	SN012423
19	SN012419	24	SN012424
20	SN012420	25	SN012425
21	SN012421	26	SN012426
		27	SN012427

# Endotracheal Tube



- Used to protect airway under general anesthesia congenital malformation and disease of the upper airway mechanical ventilation perinatal resuscitation and various form of respiratory disorder.
- Made from non-toxic, non-irritant medical grade clear P.V.C.
- Kink resistant, thin-walled, thermo-sensitive tube, softens at body temperature.
- Smooth inside / outside surface for ease in intubation, extubation and suctioning.
- Smooth polished bevel and murphy eye and tapered balloon end for atraumatic intubation and extubation.
- Universal 15mm connector at proxima end.

## TECHNICAL SPECIFICATION

### (CUFFED)

Size	Ref. no.	Size	Ref. no.
2	ET011235	6.5	ET011265
2.5	ET011235	7	ET01127
3	ET011235	7.5	ET011275
3.5	ET011235	8	ET01128
4	ET01124	8.5	ET011265
4.5	ET011245	9	ET01127
5	ET01125	9.5	ET011275
5.5	ET011255	10	ET01128
6	ET01126		

### (UNCUFFED)

Size	Ref. no.	Size	Ref. no.
3	ETC011930	6.5	ETC011965
3.5	ETC011935	7	ETC01197
4	ETC01194	7.5	ETC011975
4.5	ETC011945	8	ETC01198
5	ETC01195	8.5	ETC011985
5.5	ETC011955	9	ETC01199
6	ETC011960	9.5	ETC011995
		10	ETC011910

# Infant Mucus Extractor



- Mucus extractor is used for aspiration of secretion from oropharynx in infants to ensure free respiration.
- Funnel shaped adapter at the proximal end of suction tube to ensure easy & secured connection.
- Graduated & transparent chamber allows correct measurement & visual checking of aspirate.
- Soft, frosted and kink resistant PVC tubing.
- Spare cap is provided for safe storage of specimen for analytical testing and safe disposal.
- Sterile/ Non-Pyrogenic/ Single use.

**CAPACITY :** 25mm

**SUCTION TUBE LENGTH AVAILABLE :** 40 & 50cm

## TECHNICAL SPECIFICATION

Size	Ref. no.
8FG	IME01698
10FG	IME016910
12FG	IME016912
14FG	IME016914

# Guedel Airway



- Smooth integrated design for optimum patient comfortable.
- Safe hard non-toxic polyethylene material with color code.

## TECHNICAL SPECIFICATION

Size	Ref. no.	Size	Ref. no.
0	GA01080	3	GA01083
00	GA010800	4	GA01084
000	GA0108000	5	GA01085
1	GA01081	6	GA01086
2	GA01082		




# PVC Manual Resuscitator (AMBU BAG)



- Single patient use only.
- The set includes oxygen reservoir bag, PVC mask and oxygen tubing.
- 100% Medical level PVC material.
- Latex free components.
- All-in-one intake valve connects to oxygen reservoir bag directly.
- Additional accessories and customized label are available.

## TECHNICAL SPECIFICATION

Size	Ref. no.
ADULT	AB0120A
PAEDIATRIC	AB0120P

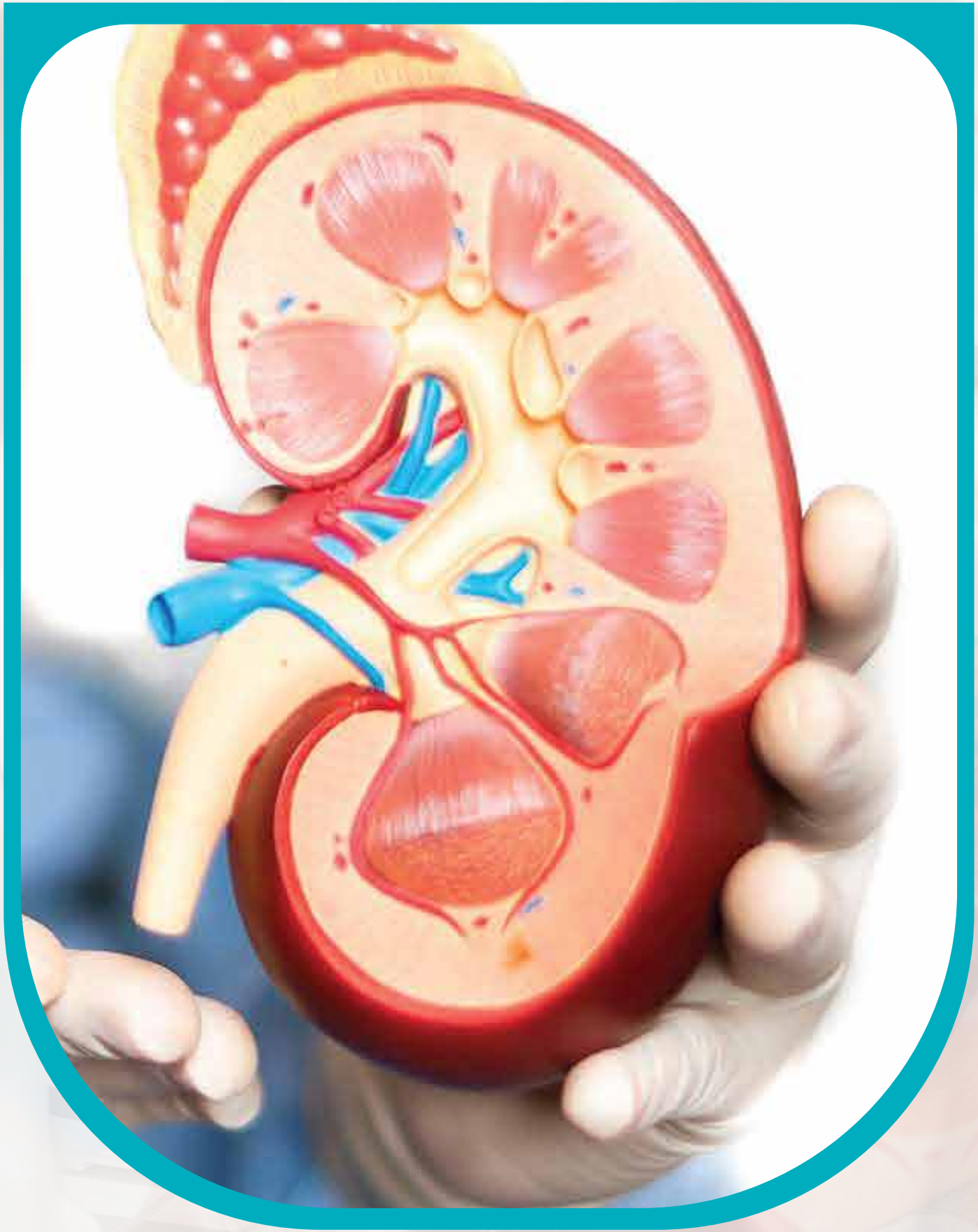

# UROLOGY

# Foley Balloon Catheter

## (Two Way)



- Made of natural latex rubber.
- Urinary tract catheter for short or long term urine drainage.
- Hard non return valve ensures trouble free inflation and deflation of balloon.
- Urine output monitoring in a critically ill or injured person.
- EO Sterile.

**AVAILABLE IN:** Paper Pouch & Box packing

### TECHNICAL SPECIFICATION

Size in FG	Balloon Capacity	Ref. no.	Size in FG	Ballon Capacity	Ref. no.
6	1.5cc-3cc	FBC01466	24	30cc-50cc	FBC014624
8	3cc-5cc	FBC01468	26	30cc-50cc	FBC014626
10	3cc-5cc	FBC014610			
12	5cc-15cc	FBC014612			
14	15cc-30cc	FBC014614			
16	30cc-50cc	FBC014616			
18	30cc-50cc	FBC014618			
20	30cc-50cc	FBC014620			
22	30cc-50cc	FBC014622			

# Three Way Foley Balloon Catheter



- Made of natural latex rubber.
- Urinary tract catheter for short or long term urine drainage.
- Hard non return valve ensures trouble free inflation and deflation of balloon.
- Urine output monitoring in a critically ill or injured person.
- EO Sterile.

**AVAILABLE IN:** Paper Pouch & Box packing

## TECHNICAL SPECIFICATION

Size	Ref. no.	Size	Ref. no.
6	FBC01216	12	FBC012112
8	FBC01218	14	FBC012114
10	FBC012110	16	FBC012116
		18	FBC012118

# Urine Collection Bag With Measured Volume Set



- Manufactured from Medical grade material.
- Designed for accurate measurement of hourly urine output during examination of the patient.
- Urine Bag is manufactured from extra strong sheeting to withstand during the long period.
- Urine bag attached directly with measured volume meter ensures auto over flow and convenience to empty the measured volume meter in the bag.
- Bag Capacity: 2000 ml. Approx.
- Chamber Capacity: 250 ml. Approx.
- Provided with "T" Type bottom Outlet is provided for easy & Quick Emptying of urine bags.
- EO Sterile.

**AVAILABLE IN:** HM Poly Pouch/Ribbon Poly Pouch.

## TECHNICAL SPECIFICATION

Size	Capacity	Packing Inner/Outer	Ref. no.
Urine Bag-Std.	2000 ml	25/250	UCMV01452000



# **SURGERY & WOUND DRAINAGE**

# Pop Bandage B.P.



- Plaster of Paris Bandage B.P. consists of leno-weave gauze fabric that is coated with a blend of the Alpha and Beta forms of Calcium Sulphate Hemi Hydrate (Plaster of Paris), together with binders and accelerators.
- Immersion time 10-12 seconds.
- Fast setting time 2.5-3.5 minutes.
- Extra creamy.
- Has temperature compensating effect.
- Low powder loss.
- Reduces chances of skin maceration.
- Manufactured as per British Pharmacopia (BP) guidelines.

## TECHNICAL SPECIFICATION

Size	Ref. No.
10cm x 2.7m	PB01621027

# Elastic Adhesive Bandage B.P.




- Best Elasticity.
- Best Quality concentrated adhesive, sticks properly and gives long lasting adhesive action.
- Dose not hamper muscular & joint flexibility and ensure smooth pressure.
- High porosity- permits air circulation.
- Easy to remove, because of thin non-adhesive edges.
- Manufactured as per British Pharmacopia (BP) guidelines.

## TECHNICAL SPECIFICATION

Size	Ref. No.
10cm x 4/6m	EB01601046



# Paraffine Gauze Dressing

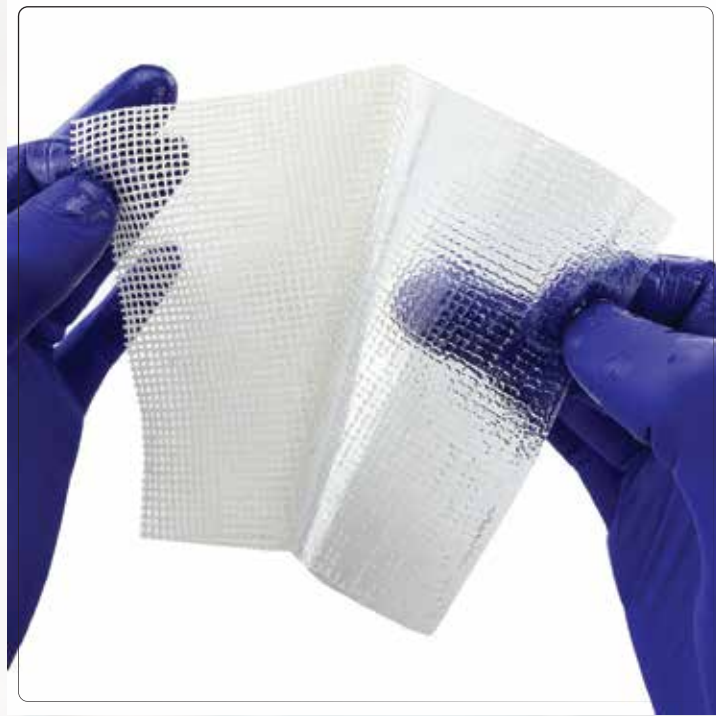

- Paraffin acts as a soothing agent. The antibiotic/antiseptic in it helps in the healing process.
- Interlocking gauze yarn threads, minimize when dressing is cut to shape.
- Does not stick to the wound and does not disturb the dressing process.
- Due to special type of butter paper loss of medicine after removal of butter paper is minimum.
- Sterilized by GAMMA radiation.

## TECHNICAL SPECIFICATION

Size	Ref. No.
10cm	PGD014910

# Chlorhexidine Gauze Dressing




- Coated with Chlorhexidine 0.5% w/w with white Soft Paraffin.
- Gauze dressing evenly impregnated with soft white paraffin containing chlorhexidine acetate B.P. 0.5%.
- Treatment For Burns, Skin Grafts, Ulcers And Other Injuries.
- Do not resterilise unused dressing or portion

## TECHNICAL SPECIFICATION

Size	Ref. No.
10cm	CGD014710

# Povidone Iodine Gauze Dressing

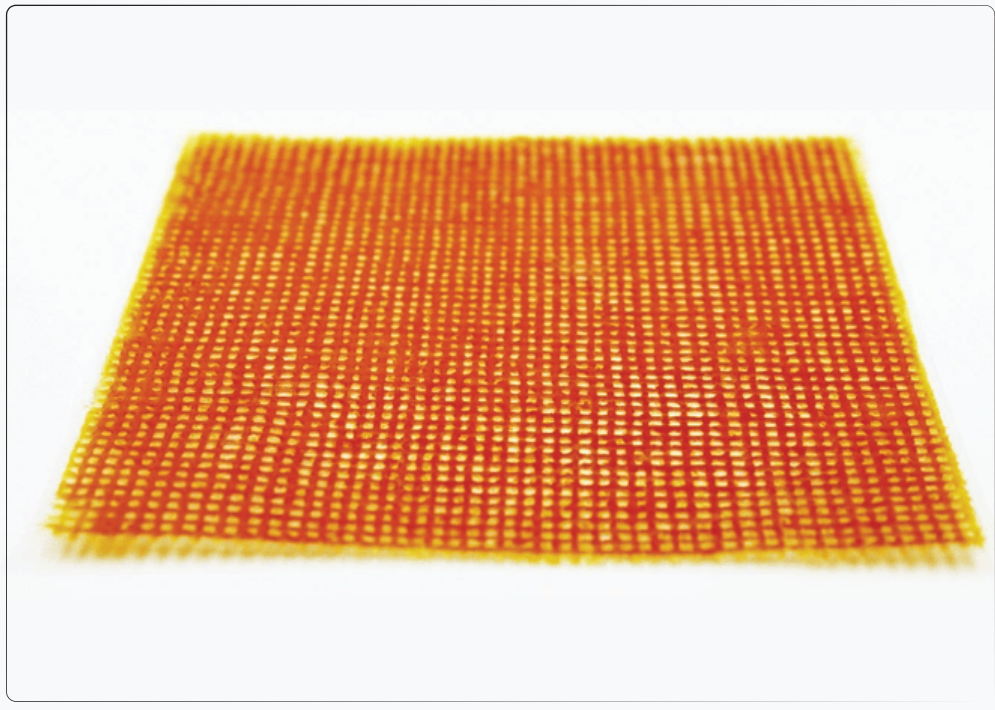

- Highly Absorbent
- Povitulle TM is an antiseptic, antimicrobial tulle dressing brownish in color, consist of knitted weave viscose fabric evenly impregnated with ointment containing.
- Soothes, Non -adherent, in presence of wound fluid POVIDONE molecule sustain releases the IODINE to come in contact with micro - organism's cell to inactivate.
- The dressing fabric is woven from superior quality of viscose yarn. Interlocking construction of threads minimizes fraying while dressing to shape.
- Interleaved for easy peel off and aseptic application

## TECHNICAL SPECIFICATION

Size	Ref. No.
10cm	PIG01291010

# Surgical Paper Tape



- Used for fixation of dressing, tubing, osteomy pouches etc on patients.
- Made from micro-porous non-woven acrylic fabric.
- Uniformly coated with hypoallergenic, non-toxic, non-irritant adhesive.
- Easy to peel off without causing discomfort or skin irritation.
- High porosity allows skin to breathe, and high resistance to water and vapour prevents the tape from peeling off.

## TECHNICAL SPECIFICATION

Size	Ref. No.
2.50cm x 9.14m	SPT0165250914

# Sutures



- Surgical sutures are medical devices used to close wounds by holding tissues together with thread. They are also known as stitches.
- Sutures can help reduce the risk of infection and support the healing process

## Silk Black Braided

- Silk suture is a non-absorbable.
- sterile, surgical suture composed of an organic protein called fibroin.
- Silk Sutures are processed to remove the natural waxes and gums.
- Braided Silk is coated with silicone and is available dyed in black with Logwood Black extract.
- Silk sutures are available in range of gauge sizes and lengths, non-needed or attached to stainless steel needles of varying types, sizes and shapes.

## Polyglactine Acid

- Polyglactin is a synthetic, absorbable, braided surgical suture made from polyglycolide and co-lactide. It's used to approximate soft tissue and in general surgery.
- 90% polyglycolide and 10% co-lactide.
- Absorbs through hydrolysis, usually within 56–70 days.
- Violet or undyed

# Combine Dressing




- The non-woven pads have soft end seals to prevent development of pressure areas.
- It also provides protection from infection and padding for major surgical wounds as it has qualities of maximum absorbency.

## TECHNICAL SPECIFICATION

Size	Ref. No.
10cm x 20cm	CD01501020

# Surgical Blades

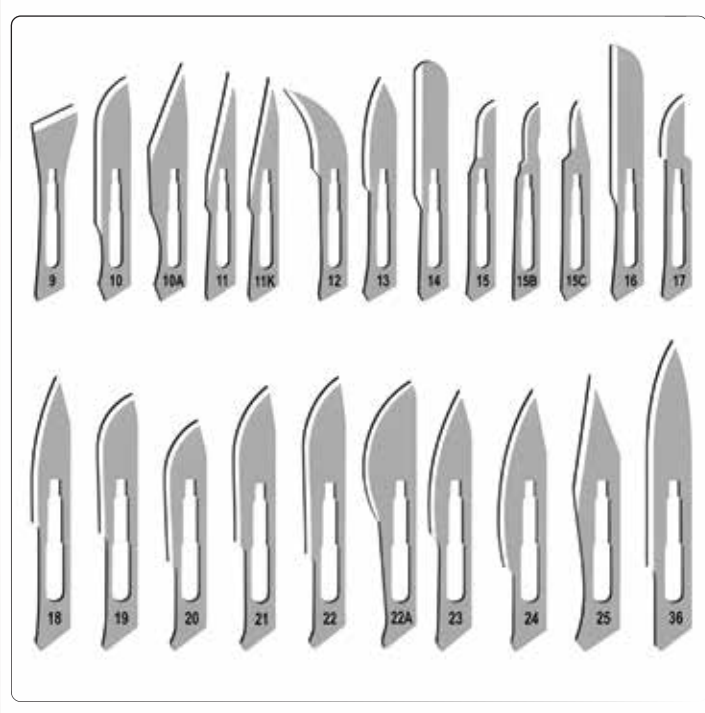

- All Popular sizes are available in both a carbon or stainless steel sterile.
- Blade shapes #9 to 17 are used in histology.
- Blade shapes #18 to 36 are used for major invasive procedures.
- Sterilization: Gamma

## TECHNICAL SPECIFICATION

Size	Ref. no.	Size	Ref. no.
10	SB012810	21	SB012821
11	SB012811	22	SB012822
12	SB012812	23	SB012823
15	SB012815	24	SB012824
20	SB012820		

# Cotton Crepe Roll



- Best Elasticity.
- Maintains uniform and controlled smooth pressure.
- Maintains Blood Circulation.
- Better aesthetic touch due to inter - woven.
- With clips.
- Reusable - Washable

## TECHNICAL SPECIFICATION

Size	Ref. No.
6 m	CB01316
8 m	CB01318
10 m	CB01310
15 m	CB01315



# Yankauer Suction Set



- Manufactured from non-toxic non-irritant medical grade PVC.
- A traumatic, soft & rounded open tip with one lateral eye.
- Provided with T type vacuum control valve with thumb control for paper maneuvering.
- Frozen surface tubing for smooth intubation.
- EO Sterile.
- Available in HM poly pack/ ribbon poly pack.

## TECHNICAL SPECIFICATION

Size	Ref. No.
2.5 m	YS011325

# Closed Wound Suction Set



- Radio opaque drainage catheter has multiple, very smooth perforations for optimum drainage, end of the perforated catheter is tip formed to minimize the trauma to patient.
- Bellow container is made of flexible plastic so that one person can easily operate to activate the section, Connecting tube is kink resistant, DEHP free with extra strength to with stand the suction.
- Bellow container is graduated marked for aspiration of 500ml/800ml fluid. Curved and siliconized alene needle is made of stainless steel for smooth sliding of the needle through tissue.
- Available in different sizes with matching size alene needle to meet moderate to heavy drainage needle, also available in mini model for small volumes after surgery (50 ml) with 6 FG & 8 FG drainage catheter & alene needle. Hanger is provided for easily handling of set. Sterile/non-pyrogenic/single

## TECHNICAL SPECIFICATION

Size	Ref. no.	Size	Ref. no.
8FG	CWS01148	16FG	CWS011416
10FG	CWS011410	18FG	CWS011418
12FG	CWS011412	20FG	CWS0114820
14FG	CWS011414		

# Underwater Seal Drainage System



- Specially designed drainage bottle suitable for pleural drainage in conjunction with intercostal drainage catheters.
- Double chamber and single chamber compact unit suitable for all cardiac and thoracic procedures.
- Range suitable for management of major, moderate and minor drainages.
- Easy to read markings help note the drained volume.
- Clearly marked initial level ensures effective under water seal.
- Write on facility on the side of scale helps in monitoring the collected drained volume precisely.
- Suction port is provided for connection with suction unit Kink resistant large bore tubing facilitates unrestricted flow.
- Metallic hangers for easy attachment to patient's bedside Sterile, individually packed.

## TECHNICAL SPECIFICATION

Size	Capacity
Adult	2000ml
Midi	1200ml
Kid	500ml



# GASTROENTEROLOGY

# Ryle Tube




- Nasogastric Feeding tube for introduction and aspiration of internal secretion.
- It has four lateral eyes to avoid blockage and easy flushing help to release any formula stuck to the inside of the tube.
- The distal end is closed by a special process to avoid any point and have marking for accurate placement.
- Soft, frosted and kink resistant P.V.C. tube radio-opaque line marked at 50, 60 & 70 cm from the tip of the accurate placement.
- Length 105 cm.
- DEHP free tube.

**Available in:** Sterile/ Non-pyrogenic/ Single use.  
With Luer connector / Without Luer connector.

## TECHNICAL SPECIFICATION

Size	Ref. No.	Size	Ref. No.
8FG	RT01168	16FG	RT011616
10FG	RT011610	18FG	RT011618
12FG	RT011612	20FG	RT011620
14FG	RT011614		

# Feeding Tube

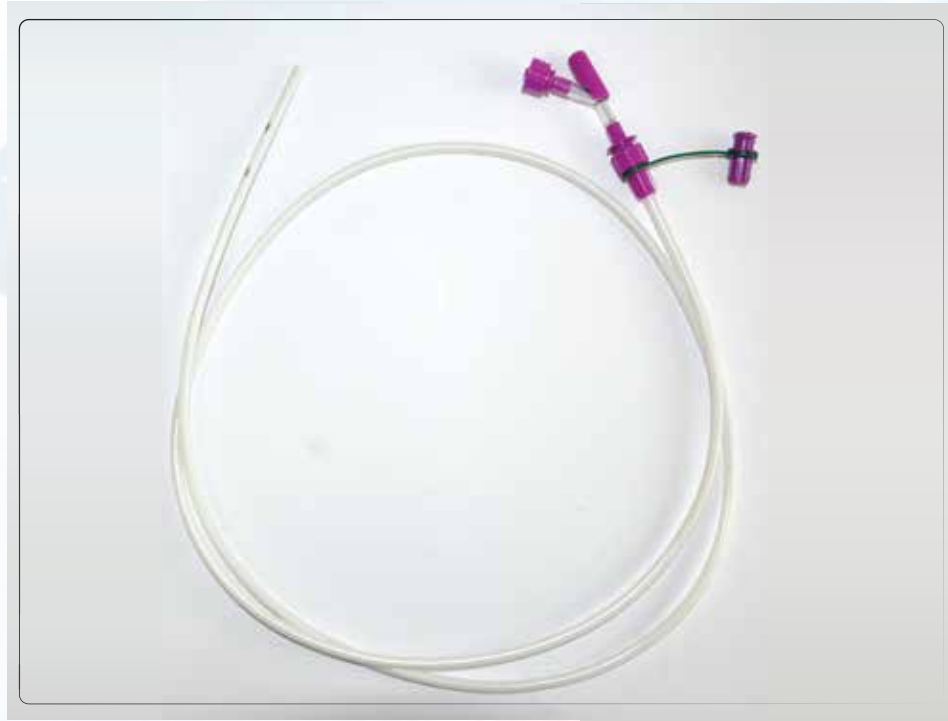

- For neonatal and paediatric nutritional feeding.
- DEHP Free tube with graduation mark after every cm.
- Atraumatic, soft rounded, closed tip with lateral eye for efficient aspiration and administration.
- X-Ray opaque line provided throughout the length of tube.
- Tube is marked at 20cm away from the distal tip.
- Closed distal end with two lateral eyes.

**Available in:** Sterile/ Non-pyrogenic/ Single use.  
 Tube length 40cm, 50cm & 100cm.

## TECHNICAL SPECIFICATION

Size	Ref. No.	Size	Ref. No.
5FG	FT01175	14FG	FT011714
6FG	FT01176	16FG	FT011716
8FG	FT01178	18FG	FT011718
12FG	FT011712		



# OTHERS

# Non-Woven Surgeon Gown



- Superior quality, non-woven fabric is obtained from reputed vendors for manufacturing these gowns.
- These gowns are duly tested on several parameters to ensure that the desired quality standards are being maintained.
- Offered gowns are both comfortable and hygienic.

**Available in:** 30,40,50,60 & 70 GSM thickness.

## TECHNICAL SPECIFICATION

Size	Ref. No.	Size	Ref. No.
Small	NWG0168S	2XL	NWG01682XL
Medium	NWG0168M	3XL	NWG01683XL
Large	NWG0168L	4XL	NWG01684XL
Extra Large	NWG0168EL		



# Non-Woven Doctors/Nurse Cap




- Adjustable soft elastic thread for comfort fittings, avoids escape of hair, dandruff and microorganisms during process.
- These caps are made from spun bond polymer non-woven fabric, spun-lace fabrics as the main raw material

## TECHNICAL SPECIFICATION

Size	Ref. No.
14"	NWC016714
16"	NWC016716
18"	NWC016718
21"	NWC016721

# Cannula Fixator

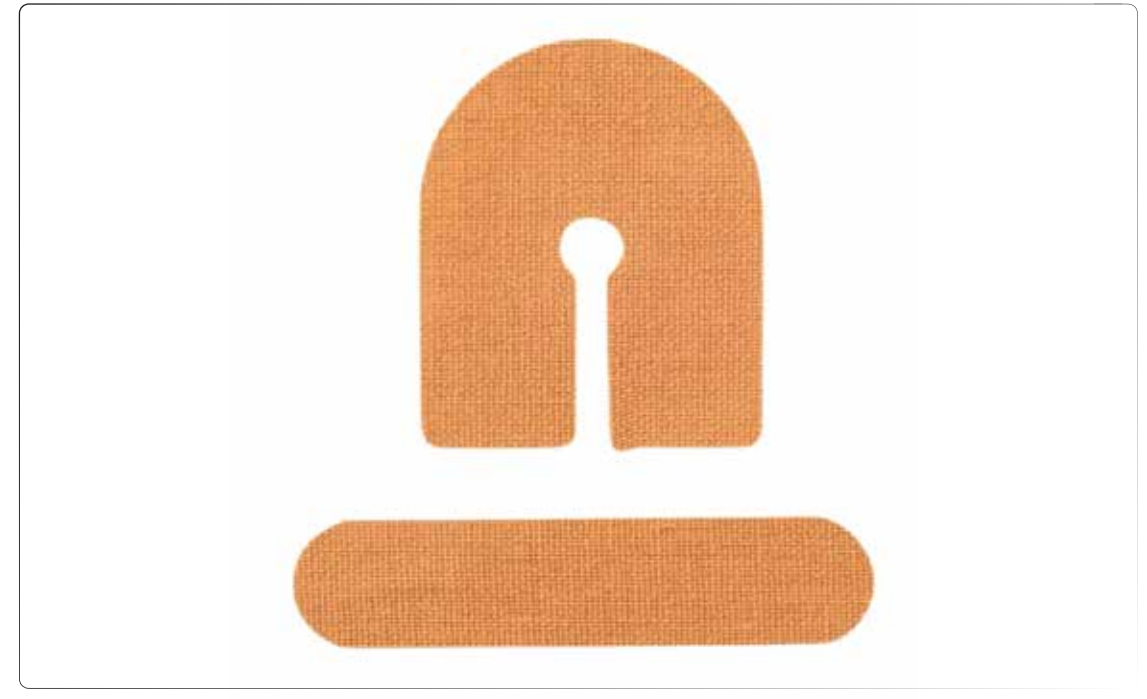

- For fixing IV Cannula
- Helps in preventing pain
- Prevents infusion related phlebitis
- Allows air circulation
- U & I Shape
- Stretches with body movements
- Manufactured as per British Pharmacopia (BP) guidelines.

## TECHNICAL SPECIFICATION

Size	Ref. No.
4.5cm x 6cm	CF0164456

# Umbilical Cord Clamp



- Designed for clamping Umbilical Cord immediately after birth of baby.
- Clamp helps stop bleeding from the blood vessels of the baby the umbilical cord.
- Provided with Non-Reopenable locks to avoid any accidental opening of clamp
- Provided with secure and leak proof jaws to ensure proper clamping.

**OPTION AVAILABLE :** Blue , White, Green

## TECHNICAL SPECIFICATION

Size	Ref. No.
5cm	UCC01705
10cm	UCC017010
15cm	UCC017015



# Gloves

# Latex Examination Gloves Non Sterile



- High Grade Natural Rubber Latex.
- Natural to White and Uniform allowing for gradual shift due to thickness variation.
- Ambidextrous, beaded cuff, smooth Storage Condition.
- The gloves shall maintain their properties when stored in a dry condition at temperature between 10° to 35° C. Avoid direct sunlight
- The gloves shall have shelf life of 5 years from the date of manufacture with the above storage condition.

**AVAILABLE :** Powdered / Powder Free

## TECHNICAL SPECIFICATION

### (POWDERED)

Size	Ref. No.
Small	LEGPP0185S
Medium	LEGPP0185M
Large	LEGPP0185L
Extra Small	LEGPP0185XS

### (POWDER FREE)

Size	Ref. No.
Small	LEGPF0186S
Medium	LEGPF0186M
Large	LEGPF0186L
Extra Small	LEGPF0186XS
Extra Large	LEGPF0186XL

# Nitrile Examination Gloves Powder Free Non Sterile



- High Grade Synthetic Latex (Acrylonitrile Butadiene Synthetic Rubber Latex).
- Natural to White and Uniform allowing for gradual shift due to thickness variation.
- The gloves shall maintain their properties when stored in a dry condition at temperature between 10° to 35° C. Avoid direct sunlight.
- Ambidextrous, beaded cuff, smooth. Storage Condition.
- The gloves shall have shelf life of 5 years from the date of manufacture with the above Storage Condition.

## TECHNICAL SPECIFICATION

Size	Ref. No.
Small	NGPF0187S
Medium	NGPF0187M
Large	NGPF0187L
Extra Small	NGPF0187XS
Extra Large	NGPF0187XL

# Sterile Latex Surgical Gloves



- High Grade Synthetic Latex (Acrylonitrile Butadiene Synthetic Rubber Latex).
- Natural to White and Uniform allowing for gradual shift due to thickness variation.
- The gloves shall maintain their properties when stored in a dry condition at temperature between 10° to 35" C. Avoid direct sunlight.
- Ambidextrous, beaded cuff, smooth. Storage Condition.
- The gloves shall have shelf life of 5 years from the date of manufacture with the above Storage Condition.

## TECHNICAL SPECIFICATION

### (POWDERED)

Size	Ref. No.
6	SGPP01886
6.5	SGPP018865
7	SGPP01887
7.5	SGPP018875
8	SGPP01888
8.5	SGPP018885

### (POWDER FREE)

Size	Ref. No.
6	SGPF01896
6.5	SGPF018965
7	SGPF01897
7.5	SGPF0189675
8	SGPF01898
8.5	SGPF018985

matu™



Excellence is Our Specialty

CE



**Matu International Pvt. Ltd.**



B-11 Jhilmil Industrial Area, Delhi -110095, India



+91-8130573898



sales@matuhealthcare.com



www.matuhealthcare.com